



FOR HIGH QUALITY MERCURY ANALYZERS

Mercury
in Liquid
Hydrocarbons

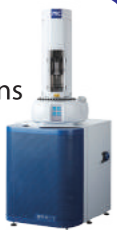

Ultra-Trace
Mercury Analyzer

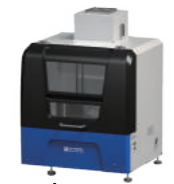

Mercury in
Water & Wastes



Fully Automated
& Continuous
Mercury Monitor



TO SERVE YOU BETTER

NIPPON INSTRUMENTS CORPORATION (NIC),
the global leader in Mercury Analyzer automation
and instrumentation, commercialized the world's first
direct thermal decomposition Mercury Analyzer back
in the 1970s. Since then, we have continued to develop,
sell, and install a vast range of Mercury Analyzers and
customized solutions around the world for regulatory
agencies and institutions, universities, research
groups, oil & gas refineries, petrochemical industries,
food providers, biology & toxicology science groups,
and more.

Thanks to our strong R&D capacities and technological
innovations, NIC is the driving force in the Mercury
Analyzer market, setting the standards for functional
scope and user-friendliness. With a wide network of
knowledgeable sales partners and qualified
service partners, we provide scientifically-sound and
environmentally responsible solutions for our customers
around the globe.

www.hg-nic.com

NIPPON INSTRUMENTS CORPORATION
14-8 Akaoji-cho, Takatsuki-shi,
569-1146 Osaka, JAPAN

Tel: +81 72 694 5195
Fax: +81 72 694 0663
E-Mail: info-nic@rigaku.co.jp



EN
MEASURING MERCURY - SINCE 1978



**HIGH-PRECISION
MERCURY ANALYZERS & MONITORS
FOR LABORATORY, RESEARCH & INDUSTRY**



YOUR WORLDWIDE PARTNER

Direct Mercury
Analyzer



Mercury in
Air & Gases



Mercury in
Flue Gases



Mercury
Survey Meters



TRUST US & ALLOW US

Version 3.0 | July 2022 | Subject to changes without prior notice.

MODEL PE-1000 Series

Highly innovative
Highly innovative, fully automated
mercury analysis for your liquid
hydrocarbons with unmatched
sensitivity and precision.
Fully compliant & accredited for
UOP Method 938



MODEL MA-3000 Series



Quick, accurate & easy to use direct mercury analyzers
for solid, liquid and gas. With no sample preparation
required and multiple configurations available to meet
your budget or automation needs.

No one does direct mercury analysis better than NIC!

MODEL AM-6F Series



Fully Automated, Continuous Mercury Monitor, the AM-6F
is designed to measure gaseous elemental mercury (GEM)
from ambient air accurately down to the sub nano-gram
per cubic-meter level. The technique used is Direct Gold
Amalgamation Sampling – CVAFS (Fluorescence
technology) Detection.

MODEL EMP-3 Series




Light-weight, compact, high sensitivity mercury survey meter,
allows direct real-time mercury measurements in the
workplace, satisfying occupational hygiene and environmental
safety guidelines like **WHO, OSHA, NIOSH, ACGIH** and more.
Real-time results are highly crucial in order to allow quick
decision making to minimize human exposure.

MODEL WA-5 Series




Ultra-trace mercury analysis in gases ranging from ambient
air to LNG / LPG by certified **EPA, ASTM, ISO & UOP**
methods. Available in AAS or AFS versions with multiple
innovative automated sampling options.

MODEL RA-4000FG+ Series

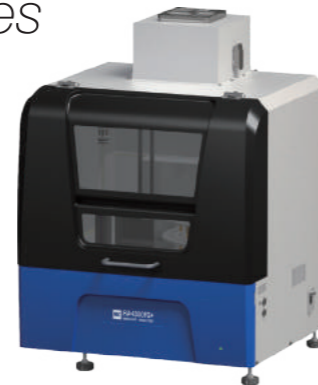

Completely automated solution for your ultra-trace mercury
(CVAFS) analysis by regulatory methods **USEPA 1631, 245.7,**
EN-1483, ISO 17852 & more. Automated addition of oxidant,
quenching agent & reduction reagent combined with
direct-purge technology in a fully enclosed “clean-room”
environment.

MODEL RA-4000A Series

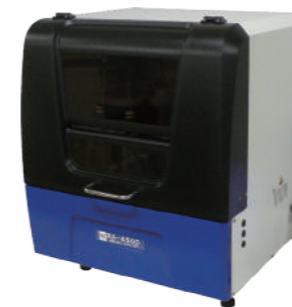

High-throughput trace level mercury (CVAAS) analyzer that
satisfies and surpasses most common regulatory permits like
USEPA 245.1, 245.2, 7470A, 7471B, APHA-311, EN-1483,
ISO-17846, ASTM D 3223, JIS K0102 and more. Available as
80-positions automated mercury analyzer or 80-positions
automated in-situ acid digestion mercury analyzer.

MODEL SGM-9 Series



In this revolutionary & innovative design, NIC has produced an
ultra-compact mercury monitor for on-site, real time
measurement of mercury in stack gases. Simple yet reliable &
lightweight for transport, with quick setup and ease of use.
A great option for real-time verification of current CEMs or
for optimization of stacks currently using absorbents for
mercury removal.